

Instruction for use – Smart Strap

Catell recommends that a licensed or otherwise qualified health professional fit and adapt the product. Users with diabetes, sensitive skin, known contact allergies, impaired sensibility or poor circulation must be extra observant and check for pinching, chafing or edema when using the brace. Always check material content below. Report any serious incidents related to the device to the manufacturer and the competent authority in the Member State where the user and/or patient resides. Contact your prescriber or Catell if questions about the product occur.

Area of use


The brace supports the user's hand and wrist after injury, surgery, chronic pain, osteoarthritis or overuse problems. The brace can also be used for preventive purposes in heavy or repetitive movements.

Application

1. Open the two hook-and-loop straps. Place your hand in the brace. **Picture 1.** Note that the picture shows right hand. For application on the left hand, rotate the brace 180° and follow the instructions below.

2. Close all hook and loop straps so that they overlap each other and the brace fits comfortably.

Picture 2, 3 and 4.

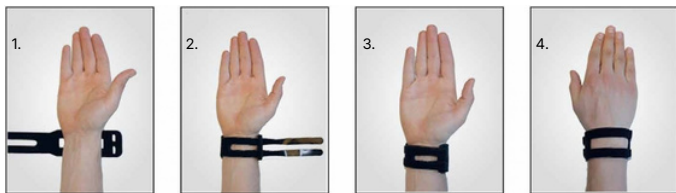

When necessary adjust the brace. It is important that the brace fits comfortably without pinching, chafing or affecting blood circulation.

Material

Neopren: synthetic rubber, polyamide. Free from latex.

Washing instructions

Close the hook and loop straps. Follow the washing instructions below. Use a washing pouch.

Let the brace drip dry.



Size

One size. The product is bilateral. Single patient use.

If the product contains neoprene, use is recommended for maximum four hours at a time.

Distributed by:

**Australian
Medical Supplies**